

**iVascular**<sup>®</sup>  
therapies for living

# nautilus

Effective large thrombus extractor (10F) to treat  
PULMONARY EMBOLISM & DEEP VEIN THROMBOSIS



Life Saving Thrombus Removal



[www.ivascular.global](http://www.ivascular.global)


# nautilus features

- > Over the wire catheter (OTW)
- > Guidewire compatibility: 0.035"
- > Coaxial luer with hemostatic valve
- > Catheter: 10F compatible
- > Usable catheter length: 120 cm
- > Proximal Shaft OD: 3.35 mm
- > Distal Shaft OD: 3.30 mm
- > Extraction area: 3.98 mm<sup>2</sup>
- > Aspiration speed: 12.00 cc/s
- > 25 cm with hydrophilic coating

## TRANSITION Catheter Design (10F)

- Rigid Proximal Shaft for Superior PUSHABILITY
- Flexible Distal Shaft for Excellent NAVIGATION
- Flat Wire Mesh tube along the entire catheter

## THE NEW STANDARD OF DELIVERABILITY

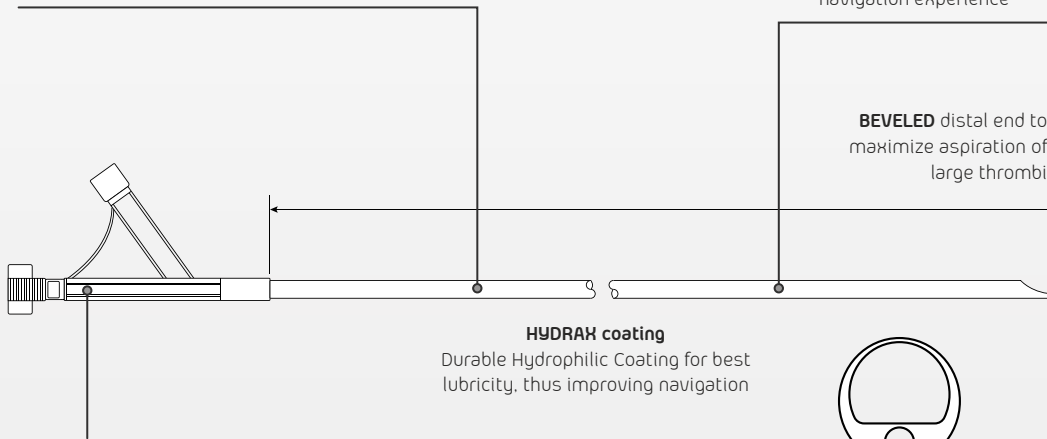

**HYDRAX coating**  
Durable Hydrophilic Coating for best lubricity, thus improving navigation

## DOUBLE LUMEN structure along the entire catheter

- Large Lumen: For Thrombus Extraction
- Small Lumen: For Guidewire
  - 0.035" GW compatibility
  - Perfect Integration with Guidewire for an EASY navigation experience

## UNIVERSAL HUB DESIGN with Hemostatic Valve

- Extraction Port: Compatible with Normal Luer Lock Syringe and Mechanical Pump
- Guidewire Port: Includes a rotating hemostatic valve to fix the wire

The biggest  
**Extraction Area** 3.98 mm<sup>2</sup>  
among competitor

## TWO metallic radiopaque markers:

- At the tip, to assess catheter position
- 15mm proximal to the tip marker, to locate the aspiration mechanism

## Ordering information

Cat No	Description
DET C35 128 100	Nautilus 10F thrombus extraction catheter
	Three-way stopcock
	Flow regulator with extension line
	60 ml vacuum syringe with sealing plunger

Distributed by:



Manufactured by:

**Life Vascular Devices Biotech, S.L.**  
www.ivascular.global  
info@ivascular.global